



# ICP Medical

**Personal Protection Apparel**  
Protection. Comfort. Infection Control.



Infection  
Control  
Products

**HELPING THOSE WHO HELP OTHERS**

## ICP Medical | Gown Styles

### Coverage Styles



Open Back



Full Coverage



Side Ties



### Neck Closure



Wrap Gown



Over the Head



Hook and Loop Closure

### Cuff Styles



Elastic Wrist w/ Thumb Loop



Thumb Loop



Knit Cuff

FOR MORE INFORMATION:

[info.ICPMedical@teammtech.com](mailto:info.ICPMedical@teammtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

## ICP Medical | Material Description

### Spunbond Polypropylene

For basic infection control, this non-woven fabric bonds fibers together to form a single layer that is appropriate only for very minimal fluid exposure.

### SMS Material

Strong and breathable, spunbond/ Meltblown/spunbond (SMS) is a multi-layer fabric composed of inner layers of meltblown polypropylene between outer layers of spunbond polypropylene that is ideal for extended wear. Light/ medium weight SMS is appropriate for low amounts of fluid and heavy weight SMS may be appropriate for moderate contact with fluids.

### Coated Polypropylene (“Poly-coated”) (PP+PE)

Soft, spunbond polypropylene is coated with a layer of polyethylene (plastic) film.

### Polyethylene (Plastic) Film (PE)

A single layer of polyethylene (plastic) film provides protection at an affordable cost.

## AAMI Gown Protection Levels

ANSI/AAMI; PB70 Barrier Performance	Test Method	Test Definition		Requirements	Anticipated fluid exposure
<b>AAMI Level 1</b>	Water resistance: Impact penetration AATCC 42”	AATCC 42	Measures the resistance of fabrics to the liquid penetration of water by impact	Water impact ≤ 4.5 g	Minimal fluid levels (some resistance to water spray)
<b>AAMI Level 2</b>	Water resistance: Impact penetration AATCC 43	AATCC 42	Measures the resistance of fabrics to the liquid penetration of water by impact	Spray Impact ≤ 1.0 g	Low fluid levels (resistant to water spray and some resistance to water penetration under constant contact with increasing pressure)
	Water resistance: Hydrostatic pressure AATCC27	AATCC 27	Measures the resistance of fabrics to the liquid penetration of water by impact under constant and increasing hydrostatic pressure	Hydrostatic Pressure ≥ 20 cm	
<b>AAMI Level 3</b>	Water resistance: Impact penetration AATCC 43	AATCC 42	Measures the resistance of fabrics to the liquid penetration of water by impact	Spray Impact ≤ 1.0 g	Moderate fluid levels (resistant to water spray and some resistance to water penetration under constant contact with increasing pressure)
	Water resistance: Hydrostatic pressure AATCC27	AATCC 27	Measures the resistance of fabrics to the liquid penetration of water by impact under constant and increasing hydrostatic pressure	Hydrostatic Pressure ≥ 50 cm	
<b>AAMI Level 4</b>	ASTM F1671, Standard Test Method for Resistance of Materials Used in Protective Clothing to Penetration by Blood-Borne Pathogens Using PHI-X174 Bacteriophage Penetration as a Test System	ASTM F1671	Measures the resistance of materials used in protective clothing to penetration by blood borne pathogens using a surrogate microbe under conditions of continuous liquid contact	Pass	Blood and viral penetration resistance (2psi)

<https://www.cdc.gov/PPEInfo/Standards/Info/ANSI/AAMIPB70Class3>

AAMI PB70: This standard is made up of 2 tests for Levels 1, 2, & 3: ASTM Spray Impact Penetration Test and Hydrostatic Head Test. For Level 4 the standard is made up of 1 test: ASTM1671, Viral Penetration of blood borne pathogens

The Spray Impact Test is performed to determine if the product is protective or nonprotective. The Spray Impact Penetration Test measures the resistance of fabrics to liquid penetration by water spray impact. Results help predict barrier performance of the isolation gown’s critical zones and indicate how well the gown will perform when fluids fall or splash onto it.

The Hydrostatic Pressure Test result indicates the level of protection from 1-3. The Hydrostatic Head Test measures the resistance of fabrics to liquid penetration by water

under constantly increasing hydrostatic pressure. Results help predict barrier performance properties of the isolation gown’s critical zones and indicate how well the gown will perform when fluid pressure is applied to it.

OSHA standard 29 CFR Part 1910.1030(d)(3)(i): Provision. When there is occupational exposure, the employer shall provide, at no cost to the employee, appropriate personal protective equipment such as, but not limited to, gloves, gowns, laboratory coats, face shields or masks and eye protection, and mouthpieces, resuscitation bags, pocket masks, or other ventilation devices. Personal protective equipment will be considered “appropriate” only if it does not permit blood or other potentially infectious materials to pass through to or reach the employee’s work clothes, street clothes, undergarments, skin, eyes, mouth, or other mucous membranes under normal conditions of use and for the duration of time which the protective equipment will be used.

FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)



## ICP Medical | AAMI Level Isolation Gowns

### AAMI Level 2 Isolation Gown: ICP-8150FC

- Easy to don over the head gown
- Thumb loop with elastic cuff offers double protection at the wrist to ensure no gap between the glove and gown
- Universal Size
- Full Coverage providing enhanced coverage clinicians can depend on to prevent contact with contaminants
- Side Ties to make tying at the waist easy
- Meets AAMI PB70:2012 Level-II standards for full coverage
- Comfortable, strong and durable SMS material
- Government contracts: GSA, DLA, MSPV
- Packaging/Quantity: 10ea/10bg/100cs



FOR MORE INFORMATION:

[info.ICPMedical@teammtech.com](mailto:info.ICPMedical@teammtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)



**ICP Medical | AAMI Level Isolation Gowns**

**Dental and AAMI Level 2 Isolation Gown: ICP-8200, ICP-8250**

- Full Coverage wrap gown provides enhanced coverage clinicians can depend on to prevent contact with contaminants
- Universal and 3XL Sizes
- Comfortable Knit Cuff
- Hook and Loop neck closure: Easily adjusts neck fit for more comfort and better coverage
- Meets AAMI PB70:2012 Level-II standards
- Comfortable, strong and durable SMS material
- Government contracts: GSA, DLA, MSPV
- Packaging/Quantity: 10ea/10bg/100cs



FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

**ICP Medical | AAMI Level 2 Isolation Gowns**

**AAMI Level 2 Isolation Gowns: ICP-8700 and ICP-8800**

- Easy to don over the head gown
- Thumb loop to ensure no gap between the glove and gown
- Universal Size
- Color: White or Yellow
- Full Coverage, providing enhanced coverage clinicians can depend on to prevent contact with contaminants
- Comfortable, strong and durable SMS material
- Meets AAMI PB70:2012 Level-II standards for full coverage
- Government contracts: DLA
- Packaging/Quantity: 10ea/10bg/100cs



FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

## ICP Medical | AAMI Level Isolation Gowns

### Poly-coated Isolation Gowns: ICP-5050YFC

- Easy to don over the head gown
- Thumb loop with elastic cuff offers double protection at the wrist to ensure no gap between the glove and gown
- Universal Size
- Full Coverage providing enhanced coverage clinicians can depend on to prevent contact with contaminants
- Constructed from Polycoated material (PPSB+PE coating) for comfort and protection
- Meets AAMI PB70:2012 Level-II standards
- Government contracts: GSA, DLA, MSPV
- Packaging/Quantity: 10ea/10bg/100cs



FOR MORE INFORMATION:

[info.ICPMedical@teammtech.com](mailto:info.ICPMedical@teammtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)



## ICP Medical | AAMI Level Isolation Gowns

### AAMI Level 3 Isolation Gown with 510k clearance: RR-300FC

- Easy to don over the head gown
- Thumb loop with elastic cuff offers double protection at the wrist to ensure no gap between the glove and gown
- Universal Size
- Full Coverage provides enhanced coverage clinicians can depend on to prevent contact with contaminants
- Meets AAMI PB70:2012 Level-III standards
- Heavy SMS material provides moderate fluid protection
- Packaging/Quantity: 10ea/10bg/100cs



FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

## ICP Medical | Open Back Gowns

### RR-200

- Easy to don over the head gown
- Thumb loop with elastic cuff offers double protection at the wrist to ensure no gap between the glove and gown
- Universal Size
- Open Back providing great coverage while keeping the wearer cool
- Breathable SMS material provides comfort and protection for low fluid settings
- Packaging/Quantity: 10ea/10bg/100cs



FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

## ICP Medical | Open Back Gowns

### ICP-8150, ICP-8300

- Easy to don over the head gown
- Thumb loop with elastic cuff offers double protection at the wrist to ensure no gap between the glove and gown
- Universal & 3XL Sizes
- Open Back providing great coverage while keeping the wearer cool
- Breathable SMS material provides comfort and protection for low fluid settings
- Government contracts: GSA, DLA, MSPV
- Packaging/Quantity: 10ea/10bg/100cs



FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)



## ICP Medical | Poly-coated Open Back Gowns

### Poly-coated Gown: ICP-5050Y

- Easy to don over the head gown
- Thumb loop with elastic cuff offers double protection at the wrist to ensure no gap between the glove and gown
- Universal Size
- Open Back providing great coverage while keeping the wearer cool
- Constructed from Polycoated material (PPSB+PE coating) for comfort and protection
- Government contracts: GSA, DLA, MSPV
- Packaging/Quantity: 10ea/10bg/100cs



FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

## ICP Medical | Poly-coated Open Back Gowns

### Poly-coated Gown: RR-300, RR-300SM, RR-300XL, RR-300XXL

- Easy to don over the head gown
- Thumb loop with elastic cuff offers double protection at the wrist to ensure no gap between the glove and gown
- Universal, Small, XL and XXL Sizes
- Open Back providing great coverage while keeping the wearer cool
- Constructed from Polycoated material
- (PPSB+PE coating) for comfort and protection
- Packaging/Quantity: 10ea/10bg/100cs



FOR MORE INFORMATION:

[info.ICPMedical@teammtech.com](mailto:info.ICPMedical@teammtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

## ICP Medical | Procedure Gowns

### ICP-7050C PP/PE Gown - Tested for Chemotherapy Drugs

- Easy to don over the head gown
- Universal Size
- Comfortable Knit Cuff: Extra soft knit cuff with a snug fit to provide maximum user comfort during extended periods of use
- Open Back coverage provides great coverage while keeping the wearer cool
- Constructed from polycoated material (PPSB+PE coating) High-quantity PE-coated material provides maximum protection against fluids
- Chemotherapy drug tested according to ASTM F739-12 with no breakthrough observed during 360 minutes of testing
  - Drug Tested
  - Doxorubicin Hydrochloride
  - Etoposide (Toposar)
  - Cyclophosphamide (Cytosan)
  - Fluorouracil
  - Paclitaxel (Taxol)
  - Thiotepa
  - Vincristine Sulfate
  - Carmustine (BCNU)
  - Cisplatin
  - Dacarbazine (DTIC)
- Government contracts: GSA, DLA, MSPV
- Packaging/Quantity: 25ea/2bg/50cs



FOR MORE INFORMATION:

[info.ICPMedical@teammtech.com](mailto:info.ICPMedical@teammtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)



## ICP Medical | Procedure Gowns

### ICP-7052FCC, ICP-7053FCC PP/PE Gowns

- Easy to don over the head gown
- XL & Universal Sizes
- Comfortable Knit Cuff: Extra soft knit cuff with a snug fit to provide maximum user comfort during extended periods of use
- Full Coverage: provides enhanced coverage clinicians can depend on to prevent contact with contaminants
- Constructed from polycoated material (PPSB+PE coating) High-quality PE-coated material provides maximum protection against fluids
- Meets AAMI PB70:2012 Level-II standards
- Government contracts: GSA, DLA, MSPV
- Packaging/Quantity: 25ea/2bg/50cs



FOR MORE INFORMATION:

[info.ICPMedical@teammtech.com](mailto:info.ICPMedical@teammtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

## ICP Medical | Procedure Gowns

### Polyethylene Gown: ICP-3845

- Easy to don over the head gown
- Universal Size
- Thumb loop at the wrist to ensure no gap between the glove and gown
- Open Back coverage provides great coverage while keeping the wearer cool
- Constructed from Polyethylene material
- Government contracts: GSA, DLA
- Packaging/Quantity: 15ea/5bx/75cs



FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

## ICP Medical | PPE Disposable Gown

### Protexer Gown Specifications

VA/DoD Item #	Description	AAMI Level	Material	Color	Size	Coverage	Cuff	Neck	Waist	Latex Free	Packaging String	Available on Gov't Contracts
ICP-8150FC	SMS Isolation Gown, Yellow, Universal, Full Coverage, AAMI 2	II	SMS	Yellow	Universal	Full Coverage	Thumbloop w/elastic wrist	Over the Head	Side Tie	Yes	10ea/10pk/100cs	GSA Advantage, MSPV, DLA
ICP-8800	SMS Isolation Gown, Yellow, Universal, Full Coverage, AAMI 2	II	SMS	Yellow	Universal	Full Coverage	Thumbloop	Over the Head	Side Tie	Yes	10ea/10pk/100cs	DLA
ICP-8700	SMS Isolation Gown, White, Full Coverage, AAMI 2	II	SMS	White	Universal	Full Coverage	Thumbloop	Over the Head	Side Tie	Yes	10ea/10pk/100cs	DLA
ICP-8150	SMS Open Back Gown, Yellow, Universal Size, Over the Head		SMS	Yellow	Universal	Open Back	Thumbloop w/elastic wrist	Over the Head	Apron Tie	Yes	10ea/10pk/100cs	GSA Advantage, MSPV, DLA
RR-200	SMS Open Back Gown, Yellow, Universal Size, Over the Head		SMS	Yellow	Universal	Open Back	Thumbloop w/elastic wrist	Over the Head	Apron Tie	Yes	10ea/10pk/100cs	
ICP-8300	MS Open Back Gown, Yellow, 3XL Size, Over the Head		SMS	Yellow	3XL	Open Back	Thumbloop w/elastic wrist	Over the Head	Apron Tie	Yes	10ea/10pk/100cs	GSA Advantage, MSPV, DLA
ICP-5050Y	Poly-coated Gown, Yellow, Universal Size, Open Back, Over the Head		Poly-coated	Yellow	Universal	Open Back	Thumbloop w/elastic wrist	Over the Head	Apron Tie	Yes	10ea/10pk/100cs	GSA Advantage, MSPV, DLA
RR-300	Poly-coated Gown, Yellow, Universal Size, Open Back, Over the Head		Poly-coated	Yellow	Universal SM XL XXL	Open Back	Thumbloop w/elastic wrist	Over the Head	Apron Tie	Yes	10ea/10pk/100cs	
ICP-5050YFC	Poly-coated Isolation Gown, Yellow, Universal, Full Coverage	II	Poly-coated	Yellow	Universal	Full Coverage	Thumbloop w/elastic wrist	Over the Head	Side Tie	Yes	10ea/10pk/100cs	GSA Advantage, MSPV, DLA
RR-300FC	SMS Isolation Gown, Yellow, Universal, Full Coverage	III	SMS	Yellow	Universal	Full Coverage	Thumbloop w/elastic wrist	Over the Head	Back Tie	Yes	10ea/10pk/100cs	
ICP-3845	Polyethylene Gown, Blue, Universal Size, Open Back		Polyethylene	Blue	Universal	Open Back	Thumbloop	Over the Head	Apron Tie	Yes	15ea/5bx/75cs	GSA Advantage, DLA

\* PB70:2012: AATCC 42-Spray Impact Penetration Test and AATCC 127 Hydrostatic Pressure Test

\* ASTM FM739-12: Chemotherapy drug test

\* Latex Free: Not made with natural rubber

FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)



## ICP Medical | PPE Disposable Gown

### Protexer Gown Specifications (continued)

VA/DoD Item #	Description	AAMI Level	Material	Color	Size	Coverage	Cuff	Neck	Waist	Latex Free	Packaging String	Available on Gov't Contracts
ICP-7050C	Procedure Gown, White, Knit Cuffs, Universal Size, Open Back, Tested for Chemo Drugs		PP/PE	White	Universal	Open Back	Knit Cuffs	Over the Head	Apron Tie	Yes	25ea/2pk/50cs	GSA Advantage, MSPV, DLA
ICP-7052FCC	Procedure Gown, White, Knit Cuffs, Extra Large, Full Coverage	II	PP/PE	White	XL	Full Coverage	Knit Cuffs	Over the Head	Back Tie	Yes	25ea/2pk/50cs	GSA Advantage, MSPV, DLA
ICP-7053FCC	Procedure Gown, White, Knit Cuffs, Universal, Full Coverage	II	PP/PE	White	Universal	Full Coverage	Knit Cuffs	Over the Head	Back Tie	Yes	25ea/2pk/50cs	GSA Advantage, MSPV, DLA
ICP-8200	Dental Gown, Yellow, Full Coverage, Velcro Neck, Knit Cuffs, Universal	II	SMS	Yellow	Universal	Full Coverage	Knit Cuffs	Hook & Loop Closure	Side Ties	Yes	10ea/10pk/100cs	GSA Advantage, MSPV, DLA
ICP-8250	Dental Gown, Yellow, Full Coverage, Velcro Neck, Knit Cuffs, 2XL	II	SMS	Yellow	3XL	Full Coverage	Knit Cuffs	Hook & Loop Closure	Side Ties	Yes	10ea/10pk/100cs	GSA Advantage, MSPV, DLA
RR-500	Chemotherapy Gown, ASTM FM739-12 drug tested, USP 800	II & USP 800	PP/PE	Purple	Universal	Full Coverage	Knit Cuffs	Hook & Loop Closure	Back Tie	Yes	25ea/2pk/50cs	

\* PB70:2012: AATCC 42-Spray Impact Penetration Test and AATCC 127 Hydrostatic Pressure Test

\* ASTM FM739-12: Chemotherapy drug test

\* Latex Free: Not made with natural rubber

FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

## ICP Medical | RR-500 Chemo Drug Tested Isolation Gown

### USP 800 Compliant /AAMI Level 2

USP 800 guidelines define stricter requirements around the type of PPE to wear when handling hazardous drugs. Chemotherapy Drug Tested Isolation Gowns are designed to comply with all USP 800 requirements and are tested against chemotherapy drugs using ASTM F739-12\*.


#### Chemotherapy Gowns

Gowns must not have seams or closures that could allow hazardous drugs to pass through. Chemotherapy Drug Tested Isolation Gown seams are tested against chemotherapy drugs using ASTM F739-12.

- Carmustine (BiCNU)\*
- Cisplatin Cyclophosphamide\*\*
- (Cytosan) Dacarbazine\*\*
- Doxorubicin Hydrochloride\*\*
- Etoposide (Toposar)\*\*
- Mitomycin C\*\*
- Paclitaxel (Taxol)\*\*
- Thio-Tepa\*\*
- Vincristine\*\*
- 5-Fluorouracil\*\*

\* 80 minutes

\*\* >480 minutes

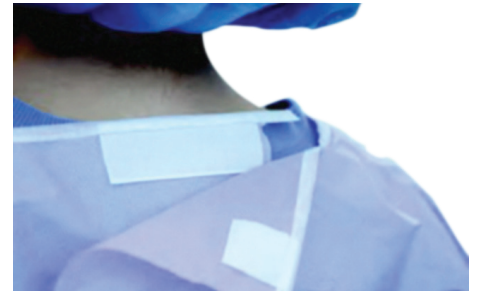

#### Gowns must close in the back

Chemotherapy Drug Tested Isolation Gowns provide fullback coverage and are equipped with long waist ties on the side of the gown to make donning as easy as possible.



#### Gowns must have closed cuffs that are elastic or knit


All gowns feature an extra soft knit cuff with a snug fit to provide maximum user comfort during extended periods of use.



#### AAMI Level 2 Fluid Protection

ICP Medical's Chemotherapy Drug Tested Isolation Gown material passes requirements for spray impact and hydrohead resistance:

- AATCC 42 Impact Penetration  $\leq$  1.0g
- AATCC 127 Hydrostatic Pressure  $\geq$  20cm
- Packaging/Quantity: 25ea/2bg/50cs

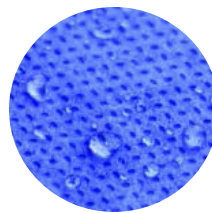

#### Unique Purple Material

Stands out from other gowns helping clinicians identify the right gown for the job without a second look. Made of PE-coated SMS material.

\* Test data available upon request.

FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

## ICP Medical | PPE Disposable Gown

### Effectively Donning and Removing Protexer Gowns

#### Helpful Hints for Donning


1. Locate label at the top of the gown. Hold at neck and allow to unfold. The patented fold allows the gown to unfold in the ideal position for donning.
2. Place head through apron style opening.
3. With thumbs facing upward, slide one arm at a time into sleeves through the thumb loop design.
4. Pull the right side of the gown behind your back and then overlap the left.
5. Gently secure ties behind the back.

Please note that any article worn on the arms, wrists or hands such as a watch or ring, could tear the gown. Please remove all adornments before donning the gown.

#### Helpful Hints for Removing

1. To remove, locate the base of the tie straps at hip and pull to release one at a time. (Do NOT pull both sides at once as you risk aerosolization of contaminants on the gown).
2. Reach up to the shoulder, grasp the gown and pull down. Repeat on the other side. (Do NOT pull both sides at once as you risk aerosolization of contaminants on the gown).
3. Roll gown into itself removing gloves into the gown being careful to only touch the inside of the gown once the gloves are removed.
4. Properly dispose of the gown.

### Removing the Gown

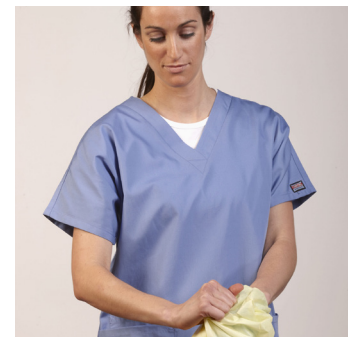

## ICP Medical | Disposable Comfort Scrubs

### Disposable Protexer Scrubs

Disposable Protexer Scrubs are a fresh alternative to laundered scrubs. The triple layer SMS breathable fabric offers supreme comfort. Whether as an alternative to laundered scrubs or as an emergency back-up supply, consider Protexer Scrubs for those inevitable occasions when theft/shrinkage, excess wear, or a surge depletes conventional inventory.

- Top: V-neck and Breast Pocket and hip pocket
- Pants: Elastic Waist
- 3-layer SMS, breathable fabric
- Color: Blue

#### Adult Protexer Scrubs

##### ICP-9000TOP

Comfort Scrub Top, Regular, Blue,  
Packaging/Quantity: 25ea/4bg/100cs

##### ICP-9000BTM

Comfort Scrub Bottom, Regular, Blue,  
Packaging/Quantity: 25ea/2bg/50cs

##### ICP-9050TOP

Comfort Scrub Top, Extra Large, Blue,  
Packaging/Quantity: 25ea/4bg/100cs

##### ICP-9050BTM

Comfort Scrub Bottom, Extra Large, Blue,  
Packaging/Quantity: 25ea/2bg/50cs

#### Pediatric Scrubs

- Top: V-neck
- Pants: Drawstring is sewn in place in the back.  
Breaks with 20 pounds of pressure.

**RR-4000S** Small, 50 sets/case

**RR-4100M** Medium, 50 sets/case



Comfort Scrubs



Pediatric Scrubs

FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)



## ICP Medical | Lab Coats

### Lab Coats

Protexer Lab Coats maximize safety and compliance for both patients and healthcare workers. The triple layer SMS material offers supreme comfort.

- 3 Pockets
- 6 Snaps
- Color: Blue
- Knit Cuff
- Traditional Collar
- Packaging/Quantity: 25ea/2bg/50cs

**ICP-9054LC** Universal

**ICP-9055LC** 2XL

- Government contracts: GSA, DLA, MSPV



FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

## ICP Medical | Bouffant Caps and Beard Covers

### Bouffant Caps

#### Description

Our bouffant caps are available in light weight single ply polypropylene or a heavy weight spunbond polypropylene. Secured with a latex free elastic band and available in multiple sizes.

Government contracts: GSA, DLA, MSPV

#### Bouffant Caps Single ply polypropylene

**ICP-CPPSB0021:** Blue, 21", 100ea/10bg/1000cs

**ICP-CP23:** Blue, 23", 100ea/5bg/500cs

**ICP-CPPSB0024:** Blue, 24", 100ea/10bg/1000cs

**ICP-CPPSB0028:** Blue, 28", 100ea/5bg/500cs

#### Heavy weight spunbond polypropylene

**ICP-CPXP24:** Blue, Extra Protection, 24", 100ea/5bx/500cs



### Beard Covers

Our beard covers are available in a lightweight single ply polypropylene material. Secured with a latex free elastic band.

**PR-BC500:** Blue, 100ea/5bg/500cs



FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)

## ICP Medical | Disposable Shoe and Boot Covers



### Shoe Covers

**ICP-SCNS0001:** Large, 100ea/3bg/300cs

**ICP-SCXT0002:** XL, 100 ea/3bg/300cs

**ICP-SCNS0003:** Universal,  
100 ea/2bg/200cs

- Elastic band for secure fit
- Fluid resistant
- Color: Blue
- Made of a polypropylene fabric and is 100% latex and silicone free
- Excellent for transitioning between "clean" and "unclean" areas
- Ideal use in medical and pharmaceutical environments, GMP facilities, pre-cleanroom gowning zones and for general hygiene use
- Anti-Skid pattern on the sole
- Non-sterile
- Government Contracts: GSA, DLA, MSPV

### PR-SCSB03

- 100ea/3bg/300cs
- Elastic band for secure fit
- Universal
- Color: Light Blue
- Ideal use for GMP or research facilities, pharmaceutical production or labs that require advanced barrier protection, especially in clean environments
- Strong fabric coated with anti-skid Polyethylene
- Made of a proprietary fabric that has higher wear resistance, fluid resistance and anti-skid properties for better foot grip
- Tested for linting compliance for cleanrooms up to class ISO 5 (FED STD 209E Class 100), very low linting
- Non-Sterile

### Boot Covers

#### ICP-BCNC0003

- Universal, 50ea/3bg/150cs
- Made of PP+PE coated fabric and is 100% latex and silicone free
- Ideal use in medical and pharmaceutical environments, GMP facilities, pre-cleanroom gowning zones and for general hygiene use
- Non-conductive strips on the sole
- Impervious
- Knee-high coverage with ankle tie
- Government contracts: GSA, DLA, MSPV

#### ICP-BCNS0002

- Universal, 50ea/3bg/150cs
- Made of PPSM+PE coated fabric and is 100% latex and silicone free
- Ideal use in medical and pharmaceutical environments, GMP facilities, pre-cleanroom gowning zones and for general hygiene use
- Traction Coating on the sole
- Impervious
- Knee-high coverage with elastic at the ankle



## ICP Medical | Stretcher Sheets

### Ultra-protective Stretcher Sheets

These disposable fluid-resistant, fitted sheets provide protection from patient secretions and body fluids. The soft-touch fabric is ultra comfortable for patients. Convenient, easy to change and sized to fit most gurneys, these fitted sheets can save valuable nursing time.

- Ultra-protective
- Minimizes risk of cross-contamination
- Fluid-repellant SMS fabric
- Elastic ends with deep corner pockets
- Tear-resistant
- Great backup for laundered sheets
- Reduces laundry/linen costs
- Latex-free

**ICP-9001GUR**

96" x 37" white SMS, 1 pc/bag;  
100 pcs/case

**ICP-9003GUR**

108" x 51" white SMS, 1 pc/bag;  
100 pcs/case

- Government contracts: GSA, DLA, MSPV



FOR MORE INFORMATION:

[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com) | 1-800-405-3044 | [www.icpmedical.com](http://www.icpmedical.com)



**TEAM Technologies** brings together a portfolio of solutions for the healthcare industry with a variety of products ranging from PPE to specialized devices. Our TEAM of experts offers end-to-end solutions for manufacturing, supply chain and distribution needs as well as products used in healthcare settings ranging from the clinic to commercial applications.

#### Our Brands/Products:

- **Rapid Refresh Privacy Curtains**
- **Protexer PPE Apparel and disposable products**
  - Gowns, Scrubs, Shoe Covers, Bouffant Caps, Lab Coats, Sleeve Covers, Beard Covers, Stretcher Sheets
- **ICP Medical PPE Organizer**
  - Customizable Configurations, Glove Boxes, Enclosed Cabinets, PPE Dispensers
- **Protexer Sticky Mats**
- **BootieButler Shoe Covers**
- **DoseLogix**
  - Topi-CLICK
- **Pure Disinfectant**
- **ICP Medical Keyboards**

#### For more information contact:

##### TEAM Technologies

5949 Commerce Blvd  
Morristown, TN 37814  
toll-free: **1-423-587-2199**  
[sales@teamtech.com](mailto:sales@teamtech.com)  
[www.teamtech.com](http://www.teamtech.com)

##### ICP Medical

13720 Rider Trail North  
St. Louis, MO 63045  
toll-free: **1-800-405-3044**  
[info.ICPMedical@teamtech.com](mailto:info.ICPMedical@teamtech.com)  
[www.icpmedical.com](http://www.icpmedical.com)

##### BootieButler

1616 Park 370 Ct.  
St. Louis MO 63042  
toll-free: **1-800-710-9863**  
[info.BootieButler@teamtech.com](mailto:info.BootieButler@teamtech.com)  
[www.bootiebutler.com](http://www.bootiebutler.com)

##### Topi-CLICK

300 Parkbrooke Place, Suite 140  
Woodstock, GA 30189  
toll-free: **1-877-870-8448**  
[info.TopiCLICK@teamtech.com](mailto:info.TopiCLICK@teamtech.com)  
[www.topi-click.com](http://www.topi-click.com)

[www.teamtech.com](http://www.teamtech.com)



Infection  
Control  
Products

